

Knee High, Anklet, Toe Piece

Measurement Form - Page 1 of 3

Patient details (Please print in black with CAPITAL letters)



Jobskin

Working together for better solutions

+61 3 9915 8000

customer.service@jobskin.com

Jobskin.com

Last name

First name

Date



Calf Sleeve - 0015



Knee Length - 0101



Anklet - 0105



Interdigital Web Spacer - 0041



Toe Piece - 1198B



Thong Toe Ankle/ Knee Length - 0539



Foot Glove - 0538

Measurements

Pleat		Pleat	
Left	Toe end	Right	
mm		mm	
	6		
	4 1/2		
	3		
	1 1/2		
	Heel 0		
	1 1/2		
	3		
	4 1/2		
	6		
	7 1/2		
	9		
	10 1/2		
	12		
	13 1/2		
	15		
Pleat	Knee End	Pleat	

Garment styles

Code	L	R	Description
0015			Calf Sleeve
0101			Knee Length
0105			Anklet (Garment length up to + 4 1/2 tape)
0041			Interdigital Web Spacer (to be worn over Foot Glove 0538)
1198B			Toe Piece ☐ Open Tips ☐ Close Tips
0539			Thong Toe Ankle ☐ Open Tips ☐ Close Tips
0539			Thong Toe Knee Length ☐ Open Tips ☐ Close Tips
0538			Foot Glove ☐ Open Tips ☐ Close Tips
0538			Individual Toe, Knee Length ☐ Open Tips ☐ Close Tips

Options and modifications

Code	L	R	Description
1169			Slant Inserts
0000			5cm Elastic
SBE			Silicon Backed Elastic ☐ 2.5 ☐ 5cm
1190			Anterior Ankle Insert (powernet & lining)
1201			Seamless Sole Standard up to the age of 10
1183			Suede Reinforcement ☐ Heel ☐ Ball
1159			Self Enclosed Toe
1160			Soft Enclosed Toe
1191			Slingback (for toe piece)
1164			Zippers (anklet)
1165			Zipper (knee-high)
ZCT			Zipper Velcro® Closure Tab (to of zipper)
1163B			Zipper Fastening Aids (along length of zip)
1172A			Zipper Loop ☐ Leather ☐ Ribbon
0000			Zipper Lining ☐ Satin ☐ Elastic
0027			Padded Zipper Lining
1187			Lining-Heel
1189			Lining-Anterior Ankle
1188			Lining-Elsewhere (specify location) ☐ Satin ☐ Cotton
1195			Fully Lined Garment ☐ Soft Satin ☐ Cotton
			SilonTEX® lining - Indicate location and size on either foot outline, TBSA map or printed photo
0027			Padding - Soft (Indicate on TBSA Map) ☐ Pocket ☐ Attached
0027			Padding - Velfoam (Indicate on TBSA Map) ☐ Pocket ☐ Attached
1202			Motif

Personalisation

Fabric Selection		Colour Selection	
☐ Conformtex		Refer to Jobskin swatch or website for fabric and thread colour choices	
☐ Powernet			
☐ Powertex			
Fabric Colour	Thread Colour	Qty	
1	1		
2	2		
3	3		

Comments

Photos assist in garment design

Leg zipper locations

	Location		Full Length	Or Specify Length (cm)
	L	R		
Lateral				
Medial				
Other				

Knee High, Anklet, Toe Piece

Measurement Form - Page 2 of 3

Patient details (Please print in black with CAPITAL letters)

Last name

First name

Date




LEFT

☐ Open Tips
☐ Closed Tips

White boxes are
lengths and grey boxes
are circumferences


TOE CIRCUMFERENCES

Toe circumference 'B'
should be 'base of toe' & 'A'
should be 'at point of
finished length' (mm), or
widest point of toe.

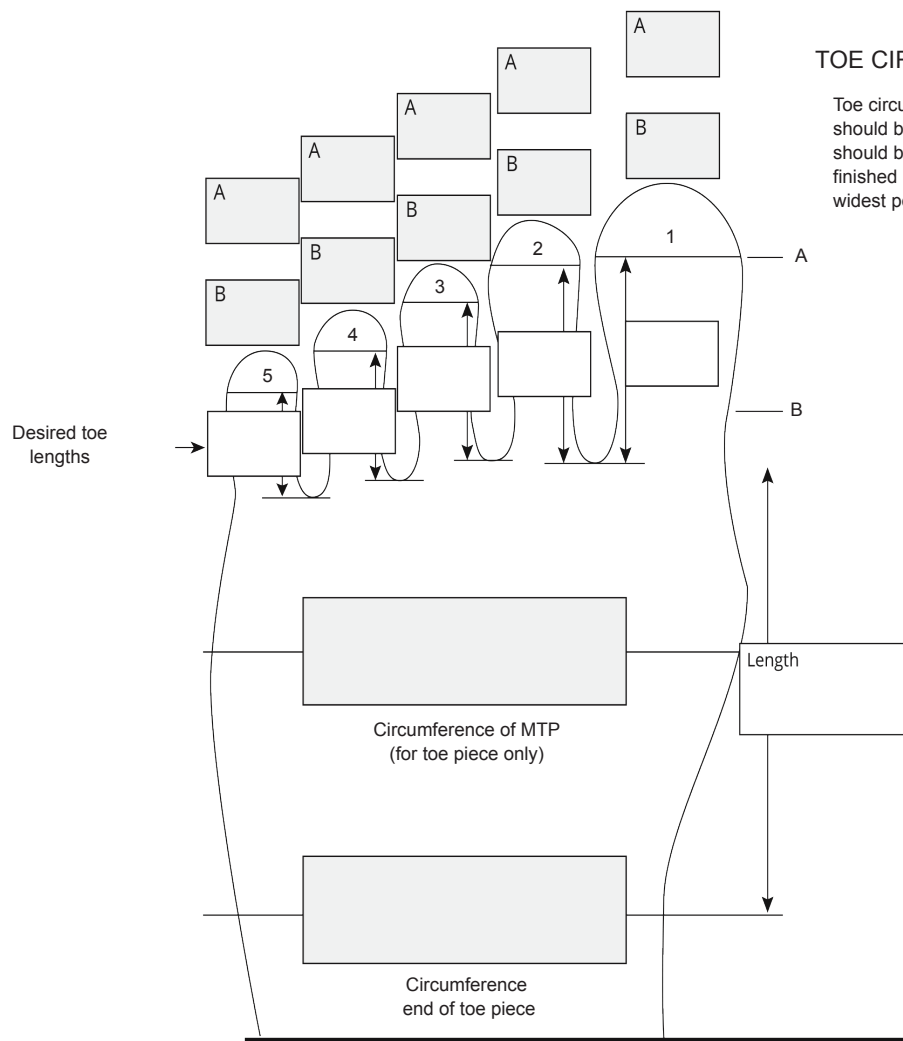

RIGHT

☐ Open Tips
☐ Closed Tips

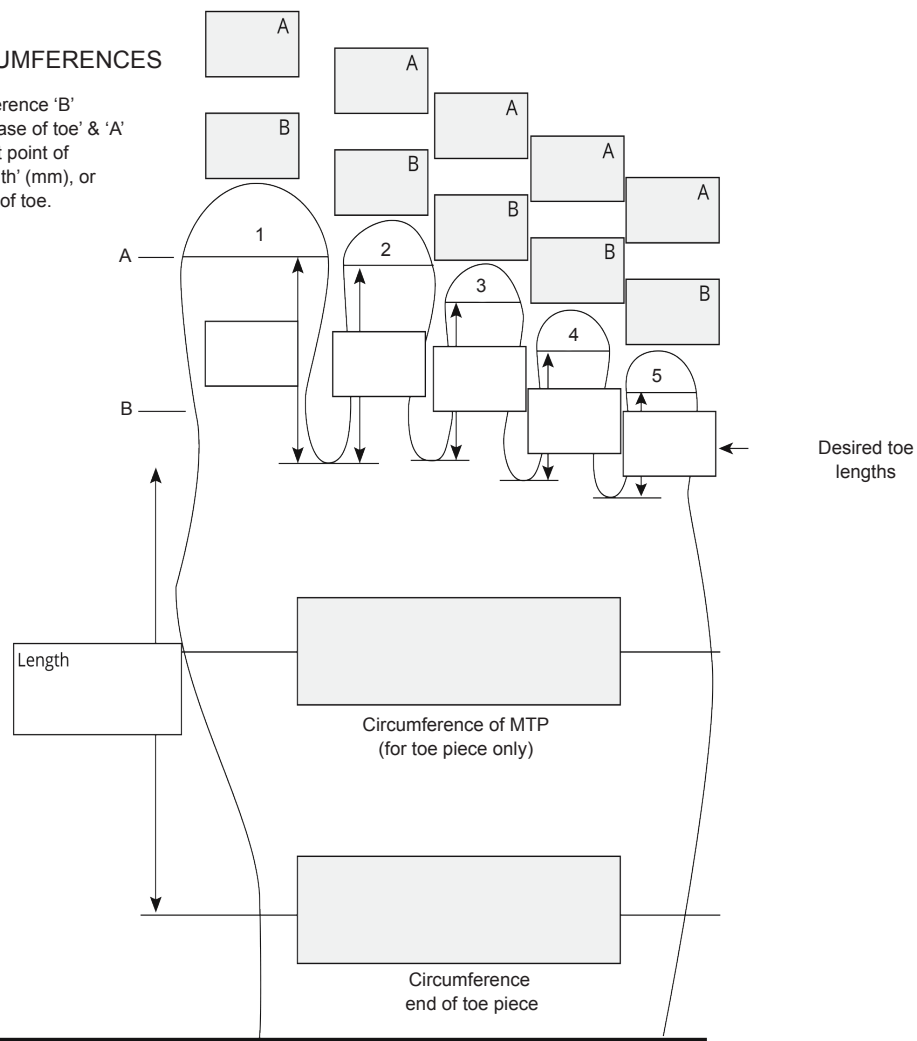

Note: FOOT OUTLINE (SEE: PG 3) MUST ACCOMPANY THIS FORM

Knee High, Anklet, Toe Piece

Measurement Form - Page 3 of 3

Patient details (Please print in black with CAPITAL letters)

Last name

First name

Date

27 26 25 24 23 22 21 20 19 18 17 16 15 14 13 12 11 10 9 8 7 6 5 4 3 2 1 0 cm

0 1 2 3 4 5 6 7 8 9 10 11 12 13 14 15 16 cm

For open tips - indicate on foot outline desired length

Please ensure all bases of webspaces are clearly marked

Place third toe on line

Toe end

Heel end

Design options




Foot outline measuring tips

Centre the patient's foot on the foot outline page as marked. Ensure that the patient is standing or sitting with the foot flat on the ground. Mark the following clearly:


- Back of the heel
- Webspaces
- Toe lengths (if measuring a toepiece or glove)
- Where the garment will finish

Use the grey JOBSKIN leg tape for these garments.
Use the pink JOBSKIN tape for paediatrics.

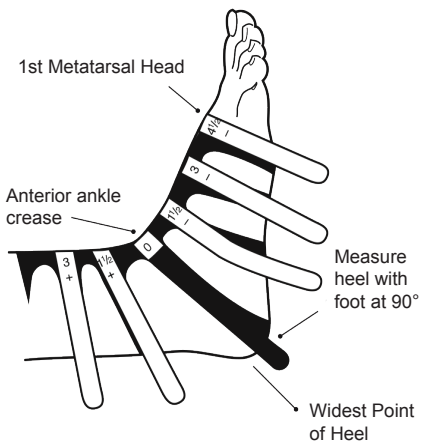


Pleat instructions

When pleating the last tape, measure the distance between the pleated tape and its neighbour from mid point to mid point (see diagram). Measurements must be recorded in millimetres only.

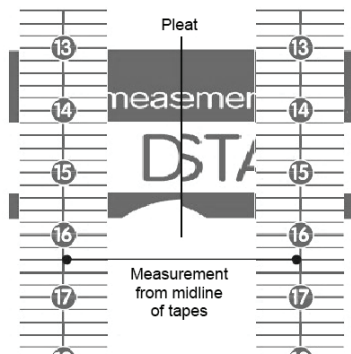

Reading the tapes

Always read along the straight edge and in direction of arrows as indicated below.

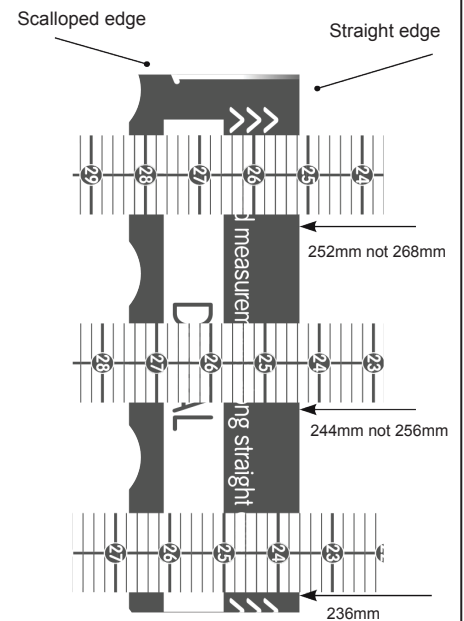

Email or return forms to:

2- 4 Queen Street
Nunawading Victoria 3131, Australia
T (Int'l): +61 3 9915 8000
E: customer.service@jobskin.com
W: www.jobskin.com

Garment returns:

2- 4 Queen Street
Nunawading Victoria 3131, Australia

Go to www.jobskin.com
for measuring videos and tips

