

LymphCrop Vest

Measurement Form - Page 1 of 2

Patient details (Please print in black with CAPITAL letters)



Jobskin

Working together for better solutions

+61 3 9915 8000

customer.service@jobskin.com

Jobskin.com

Last name

First name

Date




0561

LymphCrop Vest
Sleeveless




0562

LymphCrop Vest
2 Short Sleeves

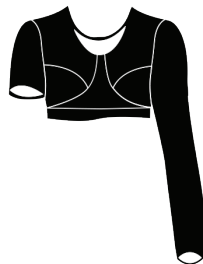

0563

LymphCrop Vest 1 Short
Sleeves, 1 Long Sleeves

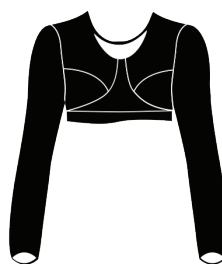


0564

LymphCrop Vest
2 Long Sleeves

The LymphCrop Vest comes standard with

- Front neckline to level of sternal component
- Back neckline with 12cm drop
- Lined front-zip opening, no tab
- 5cm silicone-backed elastic
- Lace over bra cups
- Anterior and side compression lining

Measuring points

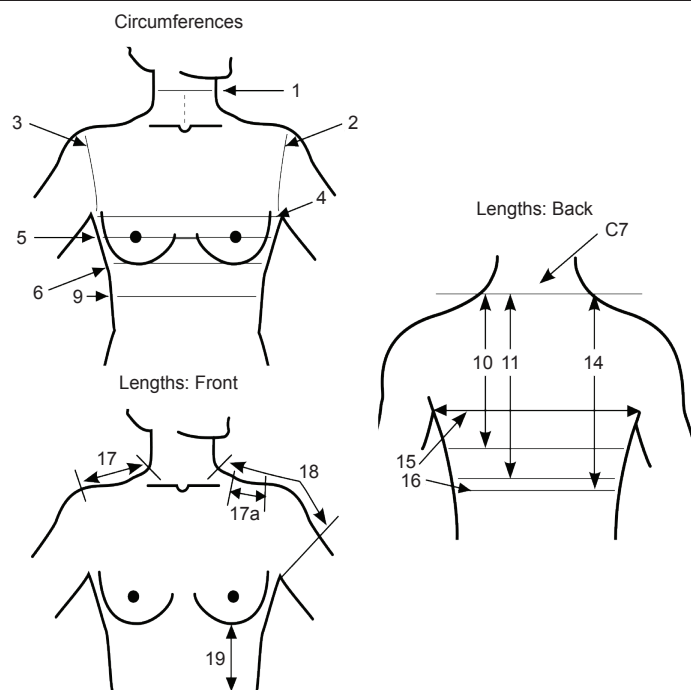

Where possible, patient should wear a bra (with prosthesis in place if required) during measurement. If not possible, position the breasts as if a bra is worn. PHOTOS TAKEN WITH PATIENT APPROVAL WILL ASSIST THE DESIGN TEAM TO PROVIDE A GOOD-FITTING GARMENT.

Garment styles

Code	Qty	Description
0561		Sleeveless LymphCrop vest
0562		Short-sleeved LymphCrop vest
0563		One short sleeve, one long sleeve LymphCrop vest
0564		Long-sleeved Lymphcrop vest

Options and modifications

Code	Qty	Description
		Prosthesis Pouch ☐ Left ☐ Right
1203		Shoulder Bracing
1190		Anterior Elbow Lining ☐ Left ☐ Right
		Replace silicone backed elastic with standard elastic (5cm)

Measurements

Neck circumference	1				
Left shoulder circumference	2				
Right shoulder circumference	3				
Chest circumference	4				
Left Chest circumference (mid sternum to mid spine)	4a				
Right Chest circumference (mid sternum to mid spine)	4b				
Circumference over nipple line	5				
Left nipple line circumference (mid sternum to mid spine)	5a				
Right nipple line circumference (mid sternum to mid spine)	5b				
Circumference immediately under breast	6				
Circumference at end of support	9				
C7 to nipple line	10				
C7 to just under breast	11				
C7 to end of support	14				
Side seam to side seam under axilla on back	15				
Side seam to side seam at end of support	16				
Shoulder width	17				
Garment shoulder width	17a				
Base of neck to top arm tape (axilla level) (Only required when sleeves are include in garment)	18				
Under breast to end of garment - front	19				

LymphCrop Vest

Measurement Form - Page 2 of 2

Patient details (Please print
in black with CAPITAL letters)



Jobskin

Working together for better solutions

+61 3 9915 8000

customer.service@jobskin.com

Jobskin.com

Last name

First name

Date

The LymphCrop Vest comes standard with

- Front neckline to level of sternal component
- Back neckline with 12cm drop
- Lined front-zip opening, no tab
- 5cm silicone-backed elastic
- Lace over bra cups
- Anterior and side compression lining
- OPTION: Contrast upper and lower bra cups rather than lace bra cups

Design choices

Neck Drop	Axilla (when sleeves includes)
☐ No Neck Drop	Powernet Axilla ☐ Left ☐ Right
	Powernet + Lining ☐ Left ☐ Right
Breast Support	Mesh Axilla ☐ Left ☐ Right
Bra Cup Size: Left Right	Satin Axilla ☐ Left ☐ Right SilonTEX® lining ☐ Indicate location and size on either TBSA map or printed photo

Upper limb circumferences

Pleat			Pleat		
Left			Right		
mm			mm		
			4 1/2		
			3		
			1 1/2		
			0		
			1 1/2		
			3		
			4 1/2		
			6		
			7 1/2		
			Elbow 9		
			10 1/2		
			12		
			13 1/2		
			15		
			16 1/2		
			18		
			19 1/2		
Pleat			Pleat		
Axilla End			Axilla End		

Only if attached armsleeve
is required

NOTE: Only pleat tapes at ends ie:at wrist and/or axilla

Personalisation

Fabric Selection		Colour Selection	
☐ Conformtex ☐ Powernet		Refer to Jobskin swatch or website for fabric and thread colour choices	
Fabric Colour	Thread Colour	Qty	
1	1		
2	2		
3	3		
Bra Cup Lace Colour			
☐ No Lace			
☐ Contrast Bra Cup			

Comments

Patient photos assist in garment design.

Design options

Contrast Bra Cup Lace



Contrast Bra Cup



Shoulder Bracing



Prosthesis Pouch



LymphCrop



Design choice definition

Alternate Neckline: The LymphCrop Vest can be made with a raised neckline. If this option is preferred, provide front and back neck drop measures to desired height.

Contrast Bra Cups: If contrast bra cups are preferred, select the décolletage colour in the colour choice box, and then select the alternate colour choice for the bra cups and the base, in the contrast colour choice box.

Shoulder bracing: A diagonal pull to help retract the shoulders to provide postural support and/or to assist with supporting the weight of large or heavily oedematous breasts.

Added Sleeves: The LymphCrop vest can be made with either short or long sleeves which enables compression for lymphoedema into the axilla. Where lymphoedema is present in both the upper limb and the chest wall, a long-sleeve may be appropriate. It is NOT recommended that a sleeve is provided for one side only, as the vest will sit asymmetrically and comfort will be compromised.

Torso measuring tips

Neck circumference: Required regardless of the height of the neckline, and measured midway between the head and the shoulders.

If neckline is to be varied from standard:

Front – measure from sternal notch point to desired neckline.

Back – measure from C7 to point of desired backline.

Bra Cup Size: Where possible, provide an accurate bra cup size based on a well-fitted bra. However, if this is not possible (particularly with challenges of fitting a bra with a lymphoedematous breast), ensure accurate 5a and 5b measures. Asking the patient to wear a tubular bandage during measuring may assist.

Shoulder circumference (2 & 3): With patient standing with hand on hip, place tape measure on top of shoulder directly above axilla, and pass through under the arm. Then lower arm down and continue tape to measure. You will notice the measure 'ease'.

Chest circumference (4): Measured at the level of the axilla patient stands with hands on hips while tape is wrapped around the circumference of the chest. Patient then lowers arms to 'ease' before reading measurement.

Nipple line (5): Circumference over nipple line should be measured with patient wearing a bra. If this is not possible, please simulate the position that the breast sits in within a bra. A tubular bandage may assist.

Side seam to side seam (15 & 16): A horizontal contoured measure from the midpoint of one axilla across the back to the midpoint of the other axilla. (15) is positioned at the height of chest circumference (4), whilst (16) is positioned at the level of the end of the garment.

If sleeve is to be included, measure upper limb with purple JOBSKIN arm tape.

If prosthesis pouch to be included, take measurements with prosthesis in place, and indicate prosthesis size.

Email or return forms to:

2- 4 Queen Street
Nunawading Victoria 3131, Australia
T (Int'l): +61 3 9915 8000
E: customer.service@jobskin.com
W: www.jobskin.com

Garment returns:

2- 4 Queen Street
Nunawading Victoria 3131, Australia

Go to www.jobskin.com
for measuring videos and tips

