



Last name	First name	Date
-----------	------------	------

Vest - 1 Short,
1 Long Sleeve **0524**Vest - 1 Sleeve,
1 Sleeveless **0526**Sleeveless Vest
0525Vest - 2 Short
Sleeves **0526**Vest - 2 Long
Sleeves **0527**Adonis 2 short
sleeves **1222**Modified Vest
0529Singlet Vest
0519

Garment styles		
Code	Qty	Description
0524		Vest with one short, one long sleeve
0526		Vest with one sleeve, one sleeveless
0525		Sleeveless vest
0526		Vest with two short sleeves
0527		Vest with two long sleeves
0523		Vest with sleeves and one gauntlet attached
0528		Vest with sleeves and two gauntlets attached
0529		Modified vest
0519		Singlet vest
0018		Torso band
1222		Adonis vest
1224		Sternal Strap - Internal
1223		Adonis - Internal

Design choices	
Neck drop	Axillas
☐ From Neck Drop cm	Powernet Axilla ☐ Left ☐ Right
☐ Back Neck Drop cm	Powernet + Lining ☐ Left ☐ Right
☐ Round Neck ☐ "V" Neck Front	Mesh Axilla ☐ Left ☐ Right
☐ Wide Neck	Satin Axilla ☐ Left ☐ Right

Nappy strap	Closures
cm	☐ Front Zipper
Waist front to waist back over nappy	☐ Back Zipper
Breast Support	☐ Velcro® Tab on zip
☐ Bra Cup ☐ Bra Cup Lace	☐ Front Velcro®
☐ Contoured Breast	☐ Back Velcro®
Cup Size	

Option and modifications		
Code	Qty	Description
1200		Funnel Neck Height: mm
0034		Turtleneck - Same Fabric ☐ Attached ☐ Detached Height: mm
0034		Turtleneck - Elastic ☐ 5cm ☐ 2.5cm ☐ 1.3cm ☐ Attached ☐ Detached
BRA CUP		Bra Cup Cup Size:
1203		Postural Support (please indicate on TBSA map)
1203		Shoulder Bracing
1190		Anterior Elbow Insert ☐ Left ☐ Right
0000		Left Elbow Contracture SeamDegrees
0000		Right Elbow Contracture SeamDegrees
1116		Irish Axilla
0001		Nappy Strap (detachable) includes Velcro tab
SBE		Silicone Backed Elastic
0027		Padding - Soft (Indicate on TBSA Map) ☐ Pocket ☐ Attached
0027		Padding - Velfoam (Indicate on TBSA Map) ☐ Pocket ☐ Attached
1163		Waistline Velcro® Attachment Tabs (Set of 4)
1202		Motif


Zippers			
Code	L	R	Description
1164			Arm Zippers
ZTC			Arm Zipper Velcro® Closure Tab
1163B			Arm Zipper Fastening Aids
0000			Arm Zipper Lining ☐ Satin ☐ Elastic
0027			Arm Padded Zipper Lining
1172A			Body Zipper Loop ☐ Leather ☐ Ribbon
1163B			Body Zipper Fastening Aids
0000			Body Zipper Lining ☐ Satin ☐ Elastic
0027			Body Zipper Padded Lining
Code	L	R	Axilla Design
			Open Mesh Standard Irish Powernet OR Powernet and Lining
			Open Mesh Standard Irish Powernet OR Powernet and Lining

Last name

First name

Date

Measuring points

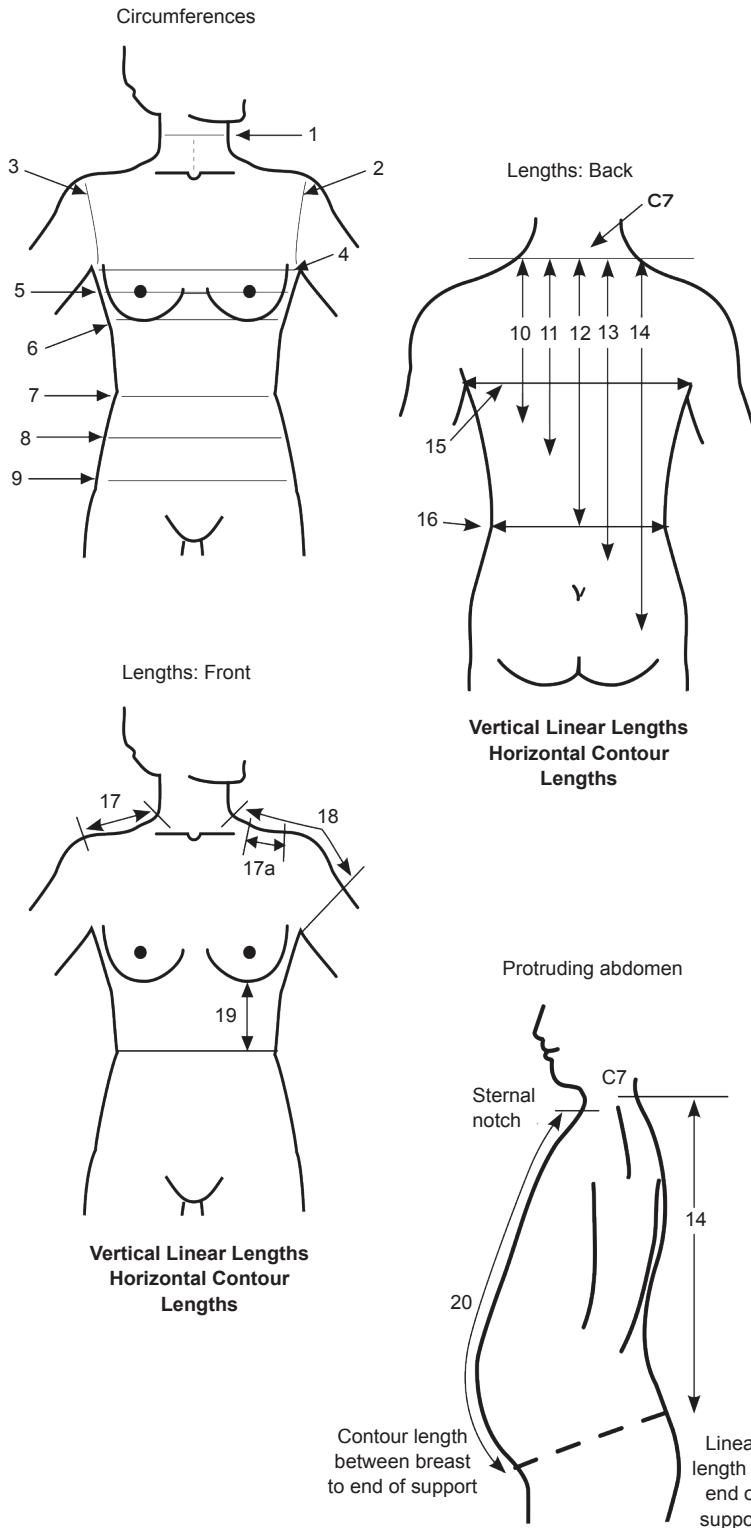

Measurements

Neck circumference	1				
Left shoulder circumference	2				
Right shoulder circumference	3				
Chest circumference	4				
Circumference over nipple line	5				
Circumference just under breast (level of breast attachment to chest)	6				
Waist circumference	7				
Circumference over largest part of belly	8				
Circumference at end of support	9				
C7 to nipple line	10				
C7 to just under breast	11				
C7 to waist	12				
C7 to largest part of belly	13				
C7 to end of support	14				
Side seam to side seam under axilla on back	15				
Side seam to side seam at waist on back	16				
Shoulder width	17				
Singlet/Sleeveless, garment shoulder width	17a				
Top arm tape to neck (arm down by side)	18				
Under breast to waist on front	19				

Protruding abdomen

Contour, sternal notch to end of support	20				
--	----	--	--	--	--

Patient to wear well fitted bra during measurement of Breast Support.

Patient details (Please print
in black with CAPITAL letters)

Last name

First name

Date

Upper limb circumferences

Pleat

Left

Pleat

Right

Wrist end

mm

mm

			4 1/2			
			3			
			1 1/2			
			0			
			1 1/2			
			3			
			4 1/2			
			6			
			7 1/2			
			Elbow 9			
			10 1/2			
			12			
			13 1/2			
			15			
			16 1/2			
			18			
			19 1/2			

Pleat

Axilla End

Pleat

Note: Only pleat tapes at ends ie:at wrist and/ or axilla

Arm zipper locations

	Location		Full Length	Or Specify Length (cm)
	L	R		
Lateral/Radial Aspect				
Medial/Ulnar Aspect				
Other:				

Personalisation

Fabric Selection

- ☐ Confortex
- ☐ Powernet

Colour Selection

 Refer to Jobskin swatch or [website](#) for
fabric and thread colour choices

Fabric Colour	Thread Colour	Qty
1	1	
2	2	
3	3	

Comments

Patient photos assist with garment design

Design options



Short-Sleeved



Long-Sleeved



Sleeveless



Adonis



Lymphcrop Vest
Use LymphCrop Vest Form



Modified Vest



Sternal Strap



Turtle Neck



Funnel Neck



Shoulder Bracing



Contoured Breast



Bra Cups (+/- lace)



Mastectomy Pouch



Irish Axilla



Mesh Axilla



Soft Axilla



Open Axilla



Nappy Strap

Design choice definition

Funnel neck: an extension of the vest up the neck without the circumferential neckline seam. Measure height as per attached powernet turtleneck.

Full neck scoop: This option will increase the width and depth of the neckline as well as decreasing compression over these areas.

Shoulder bracing: A diagonal pull to help extend the shoulders backwards and downwards to assist in improving posture as well as garment fit over the anterior chest wall.

Velcro® Closure: This closure option allows adjustability and is often used for children. It cannot be used as a front closure for women when bra cups are included.

Torso measuring tips

Neck circumference (1): Measured midway between the head and the shoulders.

Round neck back: Measure from C7 to depth of scoop.

Round neck front: Measure from sternal notch to depth of scoop.

NB. Round neck. This option is not recommended when treating the upper torso as a large scoop will significantly reduce the pressure around the shoulder, upper chest and back.

V neck front: Measure depth required from sternal notch to point of V. (Not recommended for front Velcro® opening)

Attached powernet turtleneck: Measure height from sternal notch.

Shoulder circumference (2 & 3): With patient standing with hand on hip, place tape measure on top of shoulder directly above axilla, and pass through under the arm. Then lower arm down and continue tape to shoulder. You will notice the measure 'ease'.

Arm circumference: Using Purple JOBSKIN arm tape. or Pink JOBSKIN tape for paediatrics.

For **Elbow contracture seam**, record the degrees of **deficit** in range of motion

Chest circumference (4): Measured at the level of the axilla with arms placed down to the side.

Nipple line (5): For women please provide bra cup size. Circumference over nipple line should be measured with patient wearing a bra.

Side seam to side seam (15 & 16): A horizontal contoured measure originating directly below the midpoint of one axilla across the back to the midpoint under the other axilla (15) is positioned at the height of chest circumference (4), whilst (16) is positioned at the waist.

Nappy straps: Measure over the nappy from waist at back to waist at front.

Vest length (14 & 20): We recommend the vest extend well beyond the waist for anchorage to help prevent garment riding up.

Trachea/Ostomy opening: If opening is required, it is recommended to fit the garment first, mark the position required, then return to the factory for finishing. NOTE: Ostomy openings are not recommended because of the risk of herniation. Obtain medical clearance.

Email or return forms to:

2- 4 Queen Street
Nunawading Victoria 3131, Australia
T (Int'l): +61 3 9915 8000
E: customer.service@jobskin.com
W: www.jobskin.com

Garment returns:

2- 4 Queen Street
Nunawading Victoria 3131, Australia

Go to www.jobskin.com
for measuring videos and tips

